


## #3010 Water Spray Attachment Installation Instructions

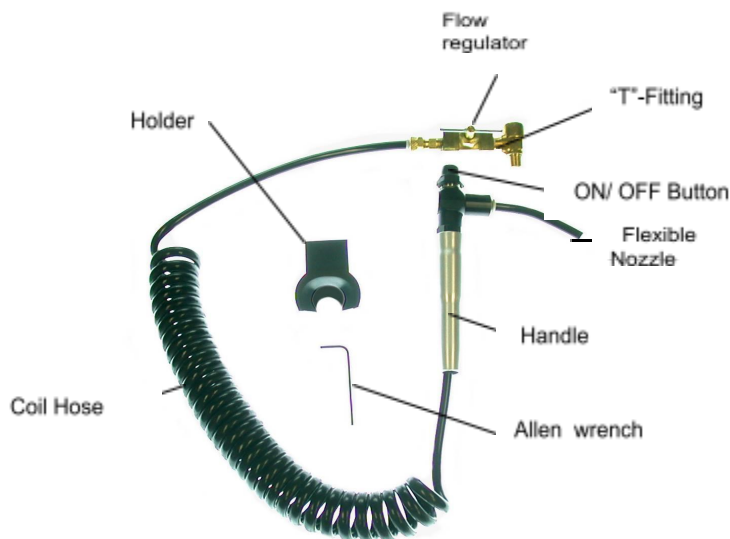

**Thanks for purchasing the #3010 WEHMER Water Spray Attachment for your Wehmer Model Trimmer.** This handy accessory is ideal for cleaning your trimming table, tools, and grinding wheels after your models are done, and for keeping your work area clean during trimming procedures.

**See the next page for installation instructions.**

## Installation Instructions

The #3010 Water Spray Attachment is mounted to the water supply line just in front of either the flow regulator or, if so equipped, the electric water control valve #3336/3337. Turn off your water supply at the source and install the “T”-fitting in the water supply line as pictured. Plumbing fittings are standard 1/8” NPT.



**Do not use Teflon thread seal tape when installing this device. If sealing of threads is required, use a pipe joint compound (pipe dope).**

The black plastic holder that stores the spray handle can be installed on the water line the flow regulator. Position the concave slot on the back of the holder over the pipe and screws from the bottom with the supplied wrench.

